


PAP Set-up Instructions

1. Remove the PAP from the packaging (plastic bag or case) and if the water chamber is in the PAP remove it by squeezing the top and pulling the chamber at the same time. If the water chamber is hard to pull out, wiggle as you pull the chamber. Once the chamber is out, fill with distilled water to the second or third line.

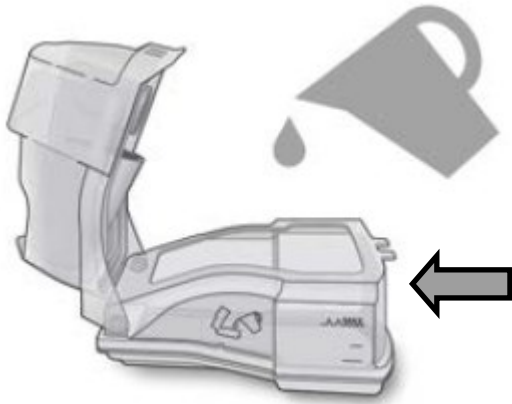

Open the water chamber by pulling out on the upper part of the chamber

2. Once the distilled water has been added, close the lid and return the water chamber to the PAP unit. Make sure you are pushing the water chamber ALL the way into the PAP unit.




DO NOT overfill the water chamber as water may enter the PAP device and tubing.


3. Hold the cuff of the air tubing and line up the air tubing connector with the connector port. Push the cuff until the connector clicks into place.




4. Plug the power connector into the rear of the PAP device. Connect one end of the power cord into the power supply unit and the other end into a power outlet.

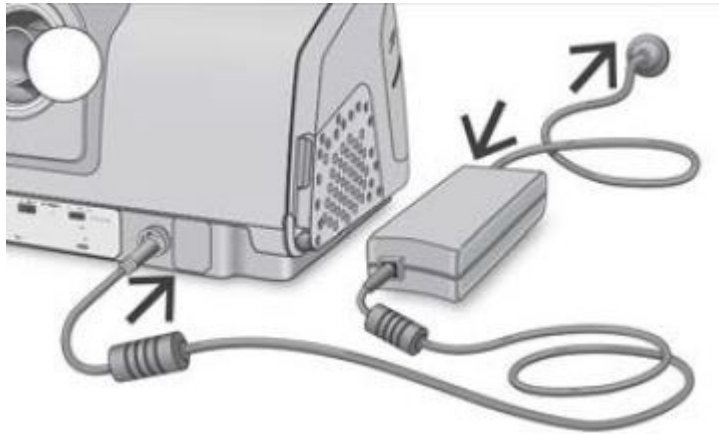

5. Connect the free end of the tubing firmly onto the assembled mask.

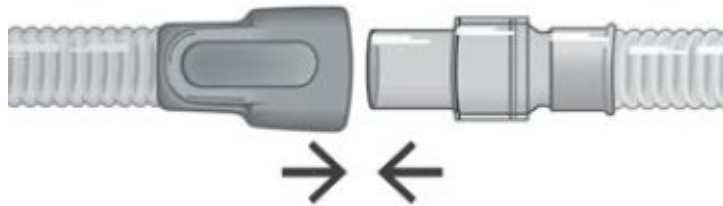

6. You are now ready to start PAP therapy by turning the unit on.

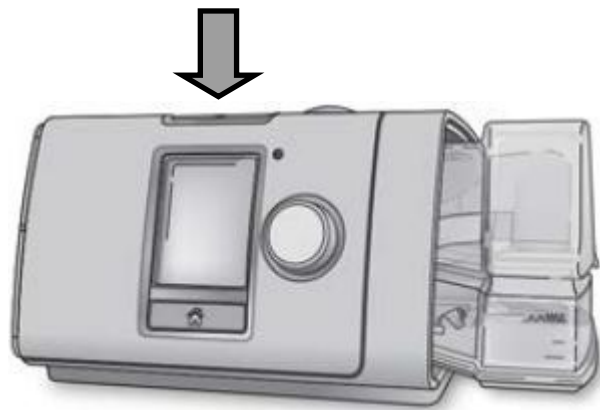

7. OFF is Green and ON is white or no lights on your PAP device.